


AMERCARE

PROTECTING YOUR PRODUCT AND ENVIRONMENT

Lead Glass Glove Port Panels

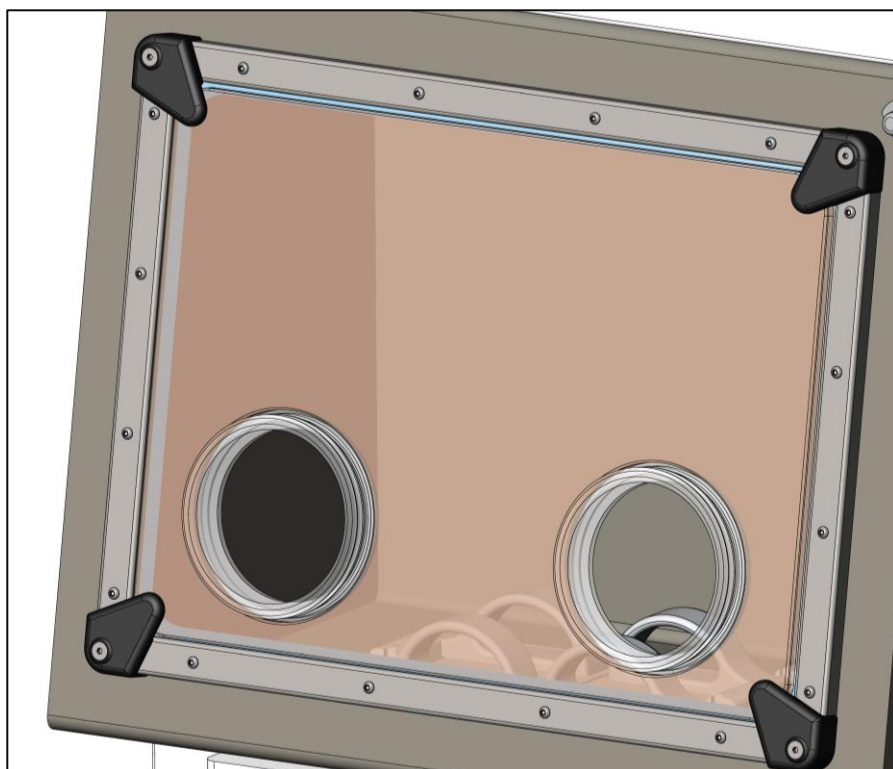

Amercure RadioPharmacy isolators are available with the option of fitting Lead Glass Glove Port panels in addition to the standard acrylic items.

These panels are available in a single thickness of nominally 8mm.


The material is transparent with a very light brown tint.

The standard acrylic panel is used to create the sealed, leak tight enclosure, the lead glass panel is then placed in front of that panel we give radiation protection to the operator and secured in place.

Specifications:

- Panels available for any isolator section
- Radiation protection to the operator – see table

Lead Glass Thickness	Approx. Lead Equivalent
Nominally 8mm	2.1mm

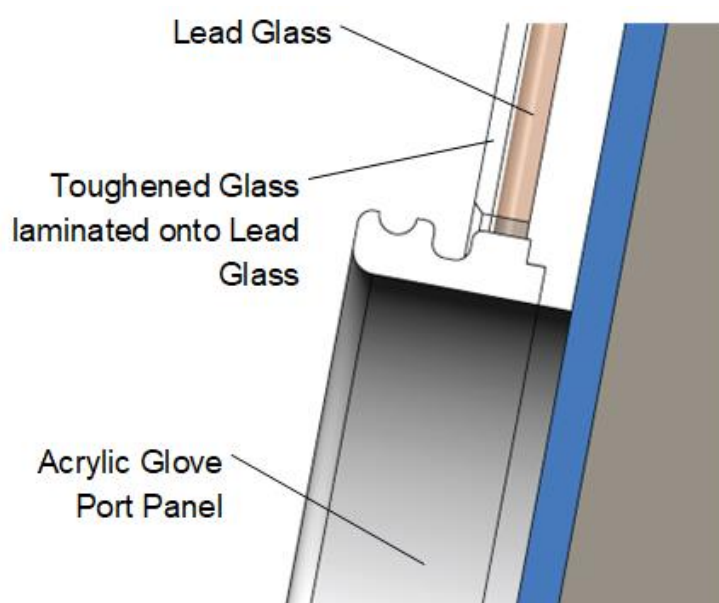

Replacement Process

Before



With just new Acrylic panel in place



With Lead Glass Panel Fitted

